



# Sensitive LC-MRM<sup>HR</sup> approach for cyclic peptide quantitation in human plasma

Lakshmanan D<sup>1</sup>, Sashank Pillai<sup>1</sup>, Rahul Baghla<sup>2</sup> and Eshani Galermo<sup>2</sup>

<sup>1</sup>SCIEX, India and <sup>2</sup>SCIEX, USA

This technical note demonstrates a sensitive method for the quantitation of pasireotide in human plasma using high-resolution accurate mass spectrometry. A lower limit of quantitation (LLOQ) of 0.05 ng/mL was achieved in extracted plasma samples (**Figure 1**).


Pasireotide is a cyclic hexapeptide therapeutic introduced in 2012 for treating Cushing's disease.<sup>1</sup> Pasireotide is a pituitary-directed medication that helps reduce cortisol secretion from the adrenal glands. Due to its key role in treating Cushing's disease, it is essential to effectively facilitate measurements of pasireotide at toxicokinetic and pharmacokinetic concentration levels in biological matrices.

This technical note demonstrates a reliable and highly sensitive workflow for supporting the quantitative analysis of pasireotide

in human plasma using a high-resolution accurate mass spectrometer.

## Key benefits for analysis of pasireotide using the ZenoTOF 7600 system

- **Sub ng/mL level of quantitation:** Achieve 0.050 ng/mL LLOQ for quantitation of pasireotide in human plasma
- **Low plasma consumption:** Reach low-level quantitation using 100 µL human plasma with increased MS/MS sampling efficiency using Zeno MRM<sup>HR</sup>
- **Effortlessly meet critical quantitative performance criteria:** Achieve accurate quantitative performance with %CV <11% at all concentration levels across a linear dynamic range (LDR) of 4.3 orders of magnitude
- **Streamlined data management:** SCIEX OS software, a 21 CFR Part 11-compliant platform, simplifies data acquisition and processing

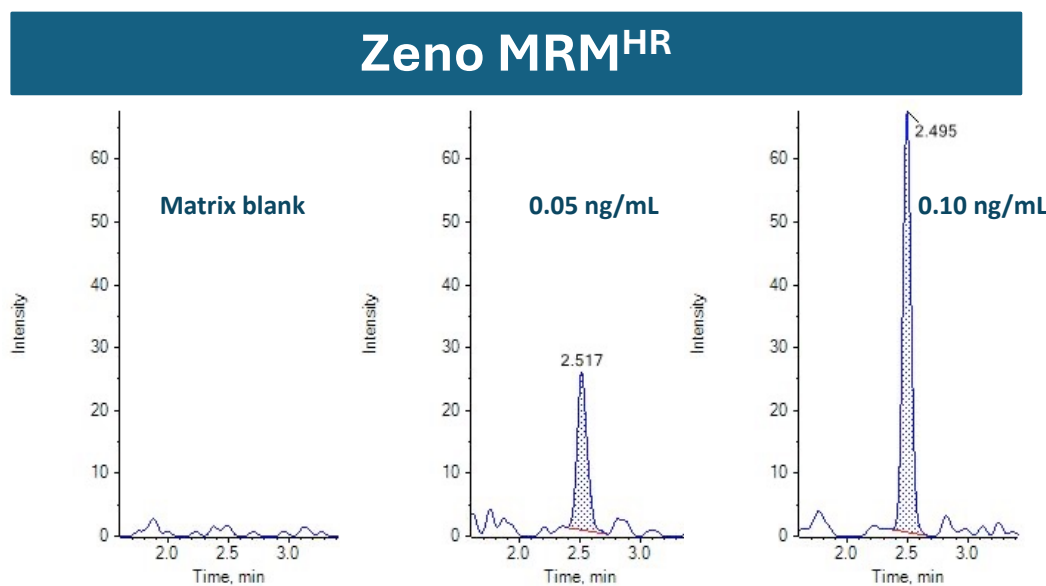

**Figure 1:** Representative extracted ion chromatograms [XICs] of matrix blank, 0.05 ng/mL (LLOQ) and 0.1 ng/mL [low quality control, LQC] for pasireotide in extracted human plasma.

## Introduction

In 2012, Pasireotide was approved for treating adult Cushing's disease patients in both the EU and USA.<sup>2</sup> It is specifically intended for those who are not eligible for surgery or have not seen improvement. Pasireotide was introduced as an alternative to somatostatin, the previous treatment option, which faced challenges due to its short half-life.<sup>3</sup> The structure of pasireotide includes a disulfide bridge, which enhances its metabolic stability and makes it suitable for a prolonged pharmacological effect.<sup>4</sup>

As a result, the complex structures of cyclic peptide therapeutics such as pasireotide require highly selective and sensitive assays to ensure accurate detection and quantitation when evaluating pharmacokinetic and pharmacodynamic effects.

## Methods

**Standard preparation:** One mg of pasireotide stock was procured from Medchem Express and dissolved in 0.1% formic acid in 50:50 [v/v] acetonitrile/water. Further, the dilutions were made in 0.1% formic acid in 20:80 [v/v] acetonitrile/water.

**Sample preparation:** Pasireotide [0.05 to 1000 ng/mL] was spiked into 100  $\mu$ L of plasma. An equal volume of 50mM TRIS buffer was added to the plasma and vortexed. The samples were subjected to solid phase extraction using a [Phenomenex Strata-X Polymeric Reverse Phase](#), 2mg 96 well plate. After loading, the sample was washed in water with 1% acetic acid in 95:5 [v/v] water/methanol.<sup>5</sup> Elution was performed with 200  $\mu$ L of methanol and dried with low nitrogen. The samples were then reconstituted with 100  $\mu$ L of 20:80 (v/v) 0.1% formic acid in acetonitrile/0.1% formic acid in water.

**Chromatography:** Analytical separation was performed on the ExionLC AE system using a [Phenomenex Kinetex C18 \[2.1  \$\times\$  100 mm, 1.7  \$\mu\$ m\]](#) column at a 0.3 mL/min flow rate. Mobile phase A was 0.1% formic acid in water and mobile phase B was 0.1% formic acid in acetonitrile. The column temperature was set to 55°C. The gradient conditions used are summarized in **Table 1**. A 20  $\mu$ L sample aliquot was injected for LC-MS/MS analysis.

| Time [min] | Mobile phase A [%] | Mobile phase B [%] |
|------------|--------------------|--------------------|
| 0.0        | 80                 | 20                 |
| 6.0        | 10                 | 90                 |
| 7.0        | 10                 | 90                 |
| 7.1        | 80                 | 20                 |
| 10.0       | 80                 | 20                 |

**Mass spectrometry:** Samples were analyzed using a ZenoTOF 7600 system operating in positive ion mode. The data were acquired using a Zeno MRM<sup>HR</sup> experiment. The optimized MS parameters are listed in **Table 2**. The summary of the Zeno MRM<sup>HR</sup> parameters is displayed in **Table 3**.

**Table 2. Source, gas and ZenoTOF 7600 system conditions.**

| Parameter              | MS       | MS/MS             |
|------------------------|----------|-------------------|
| Scan mode              | TOF MS   | MRM <sup>HR</sup> |
| Polarity               | Positive |                   |
| Gas 1                  | 50 psi   |                   |
| Gas 2                  | 70 psi   |                   |
| Curtain gas            | 35 psi   |                   |
| Source temperature     | 600°C    |                   |
| Ion spray voltage      | 5500 V   |                   |
| CAD gas                | 12       |                   |
| Declustering potential | 25 V     | See Table 3       |
| Start mass             | m/z 100  | m/z 400           |
| Stop mass              | m/z 1500 | m/z 550           |
| Q1 resolution          | NA       | Unit              |
| Accumulation time      | 0.2 s    | 0.2 s             |
| Collision energy       | 10 V     | See Table 3       |
| Zeno trap              | NA       | ON                |
| ZOD threshold          | NA       | 20,000 cps        |
| Time bins to sum       | 4        | 4                 |

**Table 1: LC gradient conditions for pasireotide.**

**Table 3: Zeno MRM<sup>HR</sup> parameters used for quantitation.**

| ID          | Precursor ion [m/z] | Fragment ion [m/z] | CE [V] | DP [V] |
|-------------|---------------------|--------------------|--------|--------|
| Pasireotide | 524.3               | 427.231            | 20     | 25     |


**Data processing:** Analysis was performed using SCIEX OS software, version 3.3.1. Peaks were integrated using the MQ4 algorithm, and a weighting of  $1/x^2$  was used for pasireotide quantitation. An XIC peak width of 0.05 Da was applied for quantitation.

### Quantitative performance on the ZenoTOF 7600 system

Cyclic peptide therapeutics such as pasireotide contain disulfide bridging, which can often introduce challenges with fragmentation. As an initial assessment, Zeno MS1 was evaluated, where precursor ion to precursor ion measurements were made without fragmentation. The approach was compared to Zeno MRM<sup>HR</sup>, where CID was applied and the product ion was measured. The matrix blank chromatograms were evaluated between Zeno MS1 and Zeno MRM<sup>HR</sup> modes. Compared to Zeno MRM<sup>HR</sup>, Zeno MS1 indicated lower selectivity with a higher background in the extracted plasma samples [Figure 2].

As a result, the Zeno MRM<sup>HR</sup> approach was used to achieve better selectivity and sensitivity for quantifying pasireotide in human plasma. The improved efficiency for total MS/MS sampling with the Zeno trap on the ZenoTOF 7600 system makes it a valuable tool for quantitative workflows that require high sensitivity and selectivity.


Using Zeno MRM<sup>HR</sup>, a calibration curve was analyzed for pasireotide at concentrations ranging from 0.05 to 1000 ng/mL. To evaluate reproducibility, each calibration standard was analyzed in triplicate. An LLOQ of 0.05 ng/mL was achieved with a minimal plasma volume of 100  $\mu$ L (Figure 1). No matrix interference was observed at the retention time of the analyte.

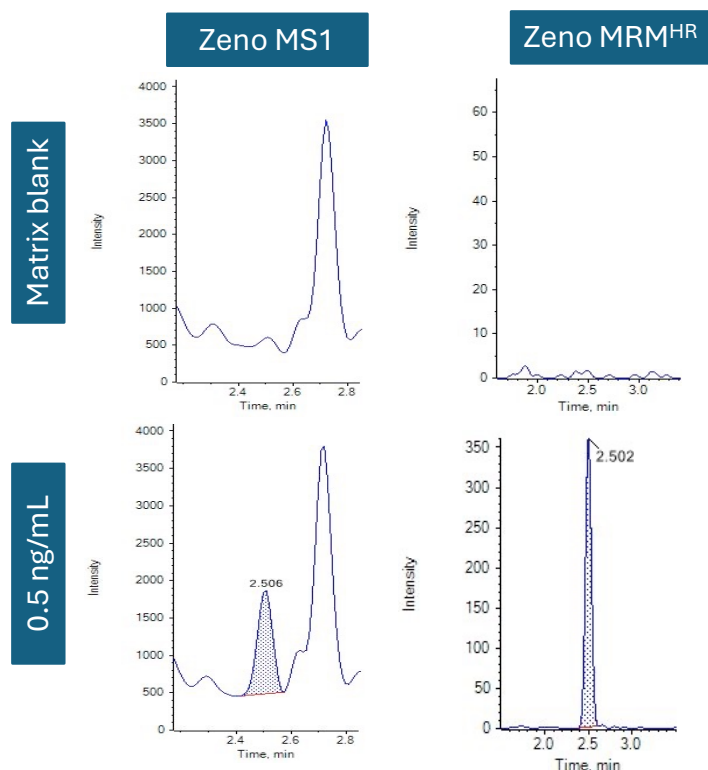

**Figure 2: Representative XICs of the matrix blank (top) and 0.5 ng/mL (bottom) of pasireotide in human plasma using Zeno MS1 (left) and Zeno MRM<sup>HR</sup> (right). A higher background was observed at the retention time of the analyte using Zeno MS1 compared to Zeno MRM<sup>HR</sup>.**

Linearity was achieved across concentrations ranging from 0.05 to 1000 ng/mL with a coefficient of determination ( $r^2$ ) of >0.994 (Figure 3), achieving an LDR of 4.3 orders of magnitude.

Calibration for PAS\_02:  $y = 3220.21070x + -6.78733$  ( $r = 0.99726$ ,  $r^2 = 0.99453$ ) (weighting:  $1/x^2$ )

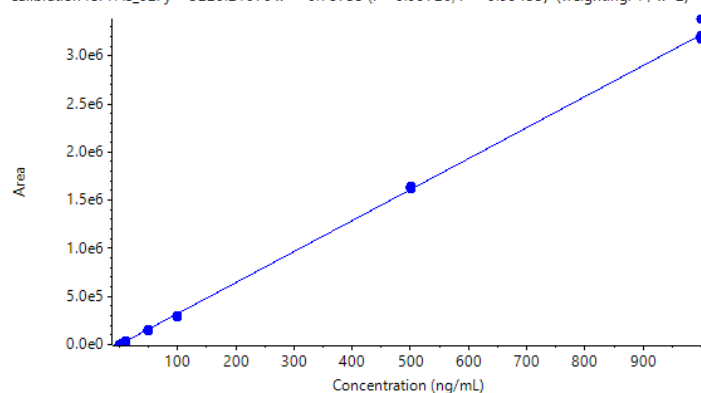

**Figure 3: Calibration curve for quantitation of pasireotide with weighing factor  $1/x^2$ .**

Analytical performance was evaluated based on the requirement that the accuracy of the calculated mean should be between 80% and 120% at the LLOQ and between 85% and 115% at higher concentrations. The %CV of the calculated mean of the concentration should be below 20% at the LLOQ and below 15% at all higher concentrations.<sup>6</sup>

For this assay, accuracy was within  $\pm 16\%$  of the nominal concentration and %CV was  $<11\%$  for pasireotide in human plasma (**Figure 4**). Calculated percent accuracy and %CV values were within the acceptance criteria at each concentration level.

| Row | Component Name | Actual Concentration | Num. Values | Mean     | Standard Deviation | Percent CV | Average Accuracy across Replicates |
|-----|----------------|----------------------|-------------|----------|--------------------|------------|------------------------------------|
| 1   | PAS_02         | 0.050                | 3 of 3      | 0.050    | 0.005              | 10.5       | 99.7                               |
| 2   | PAS_02         | 0.100                | 3 of 3      | 0.099    | 0.006              | 6.42       | 99.0                               |
| 3   | PAS_02         | 0.500                | 3 of 3      | 0.524    | 0.030              | 5.79       | 105.                               |
| 4   | PAS_02         | 1.000                | 3 of 3      | 1.062    | 0.068              | 6.44       | 106.                               |
| 5   | PAS_02         | 5.000                | 3 of 3      | 4.734    | 0.240              | 5.08       | 94.7                               |
| 6   | PAS_02         | 10.000               | 3 of 3      | 10.605   | 0.558              | 5.26       | 106.                               |
| 7   | PAS_02         | 50.000               | 3 of 3      | 47.837   | 1.465              | 3.06       | 95.7                               |
| 8   | PAS_02         | 100.000              | 3 of 3      | 91.421   | 1.488              | 1.63       | 91.4                               |
| 9   | PAS_02         | 500.000              | 3 of 3      | 506.588  | 3.180              | 0.628      | 101.                               |
| 10  | PAS_02         | 1000.000             | 3 of 3      | 1011.303 | 34.062             | 3.37       | 101.                               |

**Figure 4: Quantitative performance for pasireotide ( $m/z$  524.3  $\rightarrow$   $m/z$  427.231) analysis.** Reproducibility and accuracy results were determined from the calibration curve standards across 3 replicates at each concentration. Statistical results were summarized using the Analytics module in SCIEX OS software.

## Compliance-ready SCIEX OS software

Equivalent SCIEX OS software capabilities for regulated bioanalysis can be executed on the SCIEX 7600 system, ensuring high fidelity when performing method transfers while retaining critical compliance features.


SCIEX OS software is a closed system and requires records and signatures to be stored electronically, meeting the regulations outlined by 21 CFR Part 11. SCIEX OS software can open raw data files from any visible storage location within a closed network by using designated processing workstations. **Figure 5** illustrates the features of SCIEX OS software that are used to monitor the audit trail, acquire and process data, and configure

user access. The audit trail feature enables users to audit critical user actions and locks in data integrity. The Central Administrator Console [CAC] feature allows users to centralize acquisition and processing using a single platform to maximize efficiency for multi-instrument laboratories, independent of compliance standards. The configuration module allows users to assign roles and access as the administrator, method developer, analyst, and reviewer.



**Audit Trail**


*Easily search and filter for specific high-risk events in audit trail viewer. Built-in data integrity features allow you to tailor each functionality specifically to meet compliance needs and data security requirements.*

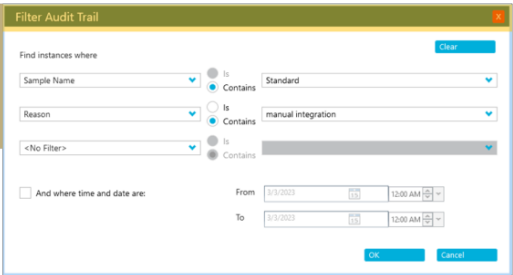



**Central Administration**


*Users can manage groups, role definitions, workstations and projects across all systems using the SCIEX OS software Central Administrator Console (CAC). It supports all regulated and non-regulated compliance standards.*

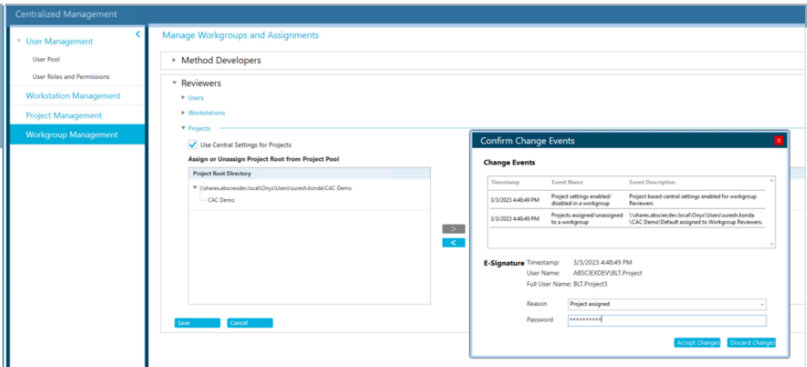



**Configuration**

*Assign users and access to administrator, method developer, analyst and reviewer roles under the audit trail module. Easily customize the role and specify level of access.*

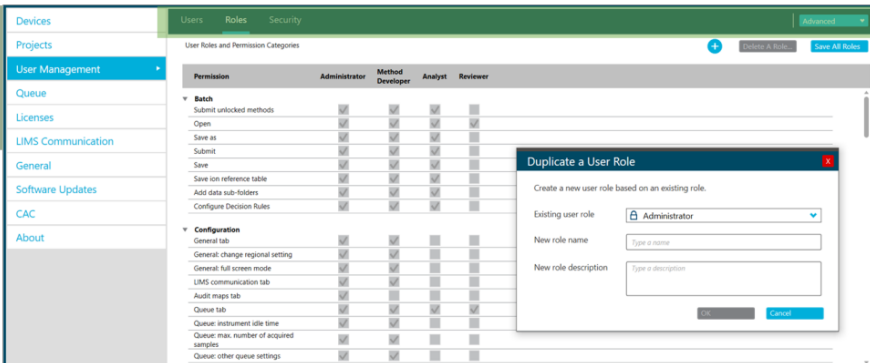

**Figure 5: Features of SCIEX OS software for monitoring user access and evaluating the audit trail.** The audit trail view allows users to filter for high-risk events easily and enables data integrity features to meet compliance requirements. The software features a Central Administrator Console [CAC] to manage users and groups, role definitions, workstations and projects across all systems. The CAC feature supports both regulated and non-regulated compliance standards. The configuration module enables users to quickly set up roles and levels of access for the administrator, method developer, analyst and reviewer levels.

## Conclusions

- An LLOQ of 0.05 ng/mL was achieved to quantify pasireotide in human plasma
- Linearity was achieved at concentrations ranging from 0.05 ng/mL to 1000 ng/mL with an  $r^2 > 0.994$  with minimal plasma volume
- Comparable quantitative performance was demonstrated with accurate and highly reproducible [%CV <11%] results on the ZenoTOF 7600 system
- A simple reverse phase solid phase extraction method was used to extract pasireotide from human plasma
- A single platform for streamlined data acquisition, processing, and management with SCIEX OS software was presented
- Retain data management and compliance-readiness [21 CFR Part 11] features using SCIEX OS software to support regulated bioanalysis on the ZenoTOF 7600 system

## References

1. Treatment of Adrenocorticotropin-Dependent Cushing's Syndrome: A Consensus Statement. [The Journal of Clinical Endocrinology & Metabolism, 2008, 2454-2462.](#)
2. [European Medicines Agency. Signifor solution for injection: summary of product characteristics. 2012.](#)
3. Pasireotide [SOM230]: Development, mechanism of action and potential applications, [Molecular and Cellular Endocrinology, 2008,286,69-74.](#)
4. A novel somatostatin mimic with broad somatotropin release inhibitory factor receptor binding and superior therapeutic potential, [Journal of Medicinal Chemistry, 2003, 46,12, 2334-2344.](#)
5. Quantitative analysis of pasireotide [SOM230], a cyclic peptide, in monkey plasma using liquid chromatography in combination with tandem mass spectrometry, [Journal of Chromatography B, 1008,242-249, 2016.](#)
6. Sensitive quantitation of insulin glargine and its metabolites in human plasma, [SCIEX technical note MKT-32026-A](#)
7. [Bioanalytical Method Validation, May 2018.](#)

The SCIEX clinical diagnostic portfolio is For In Vitro Diagnostic Use. Rx Only. Product(s) not available in all countries. For information on availability, please contact your local sales representative or refer to <https://sciex.com/diagnostics>. All other products are For Research Use Only. Not for use in Diagnostic Procedures.

Trademarks and/or registered trademarks mentioned herein, including associated logos, are the property of AB Sciex Pte. Ltd. or their respective owners in the United States and/or certain other countries [see [www.sciex.com/trademarks](http://www.sciex.com/trademarks)].

© 2024 DH Tech. Dev. Pte. Ltd. MKT-33016-A



### Headquarters

500 Old Connecticut Path | Framingham, MA 01701  
USA  
Phone 508-383-7700  
[sciex.com](https://sciex.com)

### International Sales

For our office locations please all the division headquarters or refer to our website at [sciex.com/offices](https://sciex.com/offices)