

# Sensitive and reproducible quantification of N-nitroso propranolol in a propranolol drug substance and product

Featuring a workflow for quantifying N-nitroso propranolol using the QTRAP 6500+ system

Pankaj Partani<sup>1</sup>, Sandeep Choudhary<sup>1</sup>, Preeti Bharataiya<sup>1</sup>, Upendra Gunta<sup>1</sup>, Ranjith Kumar Ponnamaneni<sup>1</sup>, Manoj Pillai<sup>1</sup>, Rahul Baghla<sup>2</sup> and Eshani Nandita<sup>2</sup>  
<sup>1</sup>SCIEX, India; <sup>2</sup>SCIEX, USA

This technical note describes the quantification of the N-nitroso propranolol impurity in a propranolol drug substance and product using the QTRAP 6500+ system. A lower limit of quantification (LLOQ) of 0.01 ng/mL was achieved for N-nitroso propranolol with high reproducibility and accuracy (Figure 1). In addition, linearity was achieved between 0.01 ng/mL and 10.00 ng/mL, providing a linear dynamic range (LDR) of 3 orders of magnitude.


Propranolol, a synthetic amino alcohol, is a competitive non-selective,  $\beta$ -adrenoreceptor antagonist extensively used to treat hypertension, angina pectoris and other cardiac diseases.<sup>1</sup>  $\beta$ -adrenergic blocking drugs such as propranolol, atenolol, metoprolol, nadolol, oxprenolol and sotalol contain amine groups. As a result, these drugs react with sodium nitrite in a hydrochloric acid solution to produce N-nitrosamines.<sup>2</sup> In March 2022, N-nitroso propranolol was reported in various

strengths of the propranolol hydrochloride extended-release capsule. Long-term exposure to N-nitroso propranolol at unsafe levels might increase cancer risk.<sup>3</sup>

Nitrosamine impurity limits are defined and must be controlled in drugs, according to ICHQ3 and M7(R1) guidelines.<sup>4-5</sup> Previous studies have revealed that the N-nitroso propranolol impurity is approximately 16-fold more carcinogenic than the N-nitroso dimethylamine (NDMA) impurity.<sup>6</sup> Therefore, there is an important need to develop assays for the sensitive detection and quantification of the N-nitroso propranolol impurity in drug substances and products.

## Key features for the quantification of the N-nitroso propranolol impurity

- **Sensitive and selective quantification:** Accomplish low-level quantification of the N-nitroso propranolol impurity using the QTRAP 6500+ system
- **Effortlessly meet critical quantitative performance criteria:** Achieve accurate and highly reproducible quantification of N-nitroso propranolol using the QTRAP 6500+ system
- **High extraction recovery:** Reach good recovery for N-nitroso propranolol in both the drug substance and drug product using a simple sample preparation workflow
- **Streamlined data management:** Data management and processing are simplified with SCIEX OS software, a 21 CFR Part 11-compliant platform

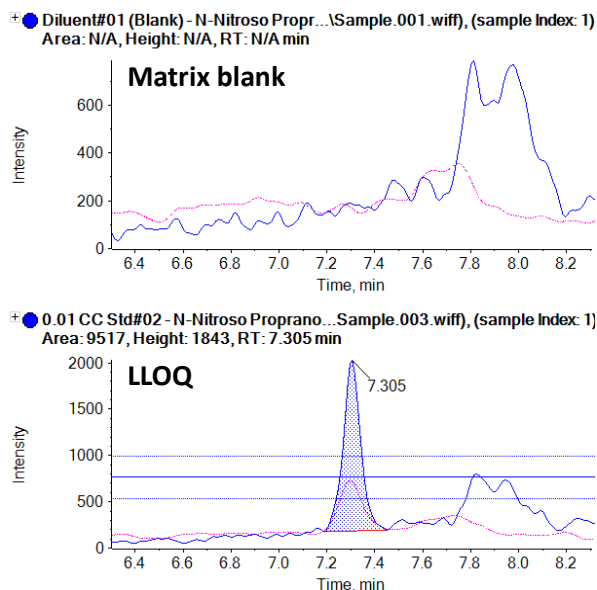

**Figure 1. Representative extracted ion chromatograms (XICs) of the quantifier and qualifier transitions for N-nitroso propranolol in the solvent blank (top) and at the LLOQ (bottom).** The LLOQ was 0.01 ng/mL. Matrix interferences were not observed at the retention time of the analyte. The ion ratio lines show the tolerance allowed of the qualifier ion to ensure accurate quantification ( $\pm 30\%$ ).

## Methods

**Standard preparation:** A 1 mg/mL N-nitroso propranolol solution was prepared in methanol. The stock solution was diluted with 80:20 (v/v), acetonitrile/water to prepare calibration standards at concentrations ranging from 0.005 to 10 ng/mL that were stored under refrigerated conditions.

**Drug substance:** Propranolol was weighed and diluted to a final concentration of 1 mg/mL in 80:20 (v/v), acetonitrile/water. The sample was shaken for 20 minutes using a mechanical wrist action shaker. After extraction, the sample was centrifuged for 10 minutes at 4000 rpm. The supernatant was collected and filtered using a 0.22 µm PVDF syringe filter. The filtered solution was transferred into an HPLC vial for analysis.

**Drug product:** The drug product was crushed and weighed. A 1 mg/mL solution of the drug product with respect to the active pharmaceutical ingredient (API) was prepared using 80:20 acetonitrile/water. The solution was vortexed for 2 minutes and shaken for 30 minutes using a mechanical wrist action shaker. After extraction, the sample was centrifuged for 10 minutes at 4000 rpm. The supernatant was collected and filtered using a 0.22 µm PVDF syringe filter. The filtered solution was transferred into an HPLC vial for analysis.

**Chromatography:** An ExionLC system was used with a [Phenomenex Kinetex biphenyl column \(3.0 x 150 mm, 2.6 µm, 100 Å\)](#) for chromatographic separation at a flow rate of 0.4 mL/min. The operating column temperature was 40°C. Mobile phase A was 1mM ammonium formate with 0.1% (v/v) formic acid in water. Mobile phase B was 0.1% (v/v) formic acid in acetonitrile. Table 1 summarizes the gradient conditions used. An injection volume of 15 µL was used for analysis.

Table 1. LC gradient.

Time (min)	Mobile phase A (%)	Mobile phase B (%)
0.0	65	35
3.0	65	35
5.0	45	55
6.0	25	75
8.0	25	75
8.1	5	95
10.0	5	95
10.1	65	35
14.0	65	35

**Mass spectrometry:** A QTRAP 6500+ system operated in positive electrospray ionization (ESI) mode was used for analysis. Table 2 summarizes the source and gas parameters used.

Table 2. Optimized source and gas parameters.

Parameter	Value
Curtain gas	40 psi
CAD gas	10
Ion spray voltage	5500 V
Temperature	450°C
Ion source gas 1	60 psi
Ion source gas 2	60 psi

The MRM conditions used for N-nitroso propranolol analysis are described in Table 3.

Table 3. MRM conditions for N-nitroso propranolol analysis.

Compound	Q1 (m/z)	Q3 (m/z)	Dwell (sec)	DP (V)	EP (V)	CE (V)	CXP (V)
N-nitroso propranolol-1	289.1	259.1	150	50	10	8	14
N-nitroso propranolol-2	289.1	72.1	150	50	10	15	15


**Note:** N-nitroso propranolol-1 was used as a quantifier transition and N-nitroso propranolol-2 was used as a qualifier transition.

**Data processing:** The data were acquired using Analyst software, version 1.7.2 and data processing was performed using SCIEX OS software, version 3.1.

## Results

ESI in positive polarity offered highly sensitive signal intensities for N-nitroso propranolol. Several chromatographic gradients and columns of various chemistries were evaluated (data not shown) to optimize the selectivity and specificity for quantifying the N-nitroso propranolol impurity in the presence of excipients, drug-related impurities and N-formyl propranolol impurities. The Phenomenex Kinetex Biphenyl LC column (3.0 x 150 mm, 2.6  $\mu$ m, 100 Å) was selected for separation, as it provided high separation power.

Under the optimal LC-MS/MS conditions, the N-nitroso propranolol impurity eluted at a retention time of 7.30 minutes. Matrix interferences were not observed in the diluent or control samples, demonstrating the selectivity and specificity of the assay. A limit of detection (LOD) of 0.005 ng/mL and LLOQ of 0.010 ng/mL for the N-nitroso propranolol impurity were achieved. The representative XICs for the matrix blank and the LOD, LLOQ and specification levels of the sample are shown

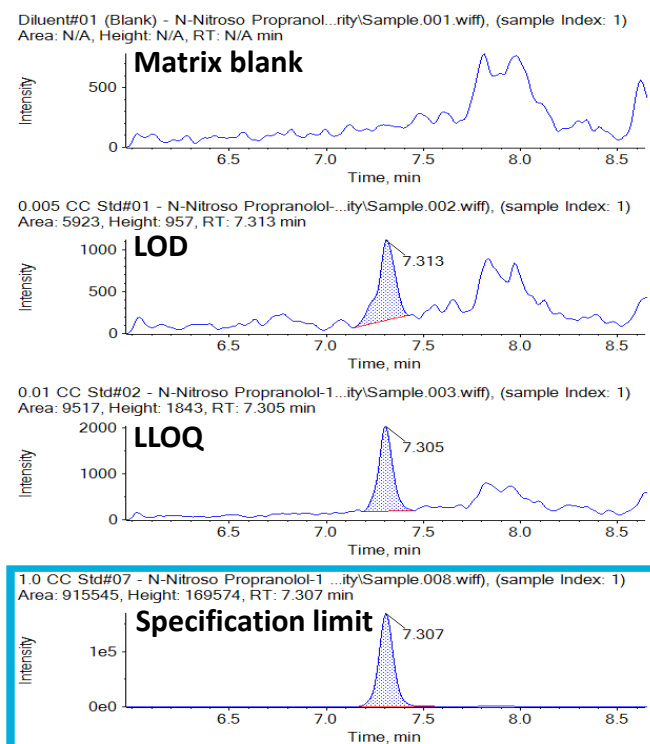

**Figure 2.** XICs of the matrix blank and sample at the LOD, LLOQ and specification limit. The LOD was 0.005 ng/mL and the LLOQ was 0.01 ng/mL. The specification limit was 1 ng/mL.


in Figure 2. The signal-to-noise (S/N) ratios observed at the LOD and LLOQ levels were >10 and >20, respectively.

The upper limit of quantification (ULOQ) was selected based on the highest levels of analytes expected to be observed and is not a reflection of system capability. The calibration curve generated for the concentrations tested is shown in Figure 3 with an overlay of the linearities for both the quantifier and qualifier transition responses. Linearity was demonstrated over a concentration range of 0.01 ng/mL to 10.00 ng/mL. A correlation coefficient ( $r^2$ ) of >0.99 was achieved for responses from both transitions.

Accuracy was expressed as percentage deviation from the nominal value at each respective concentration level. The precision of the assay was measured by the percent coefficient of variation (%CV) at each of the concentrations. Accuracy and %CV were determined by analyzing 6 replicates of the samples at the LLOQ and specification limit in a single analytical run. As shown in Table 4, the %CV values at the LLOQ and specification limit were 6.74% and 3.59%, respectively. Accuracy ranged from 96% to 115% at the LLOQ and 107% to 116% at the specification limit, averaging 105% and 110%, respectively. Overall, these data confirm the accuracy and high reproducibility of the method for quantifying N-nitroso propranolol (Table 4).

**Table 4.** Accuracy and %CV at 3 concentration levels.

Replicate No.	LLOQ (0.01 ng/mL)	0.03 ng/mL	Specification limit (1 ng/mL)
1	0.011	0.029	1.156
2	0.010	0.031	1.096
3	0.010	0.030	1.084
4	0.011	0.031	1.070
5	0.011	0.032	1.082
6	0.010	0.031	1.160
Mean	0.011	0.031	1.108
SD (%)	0.001	0.001	0.040
%CV	6.74	3.30	3.59
% Accuracy	105	103	111

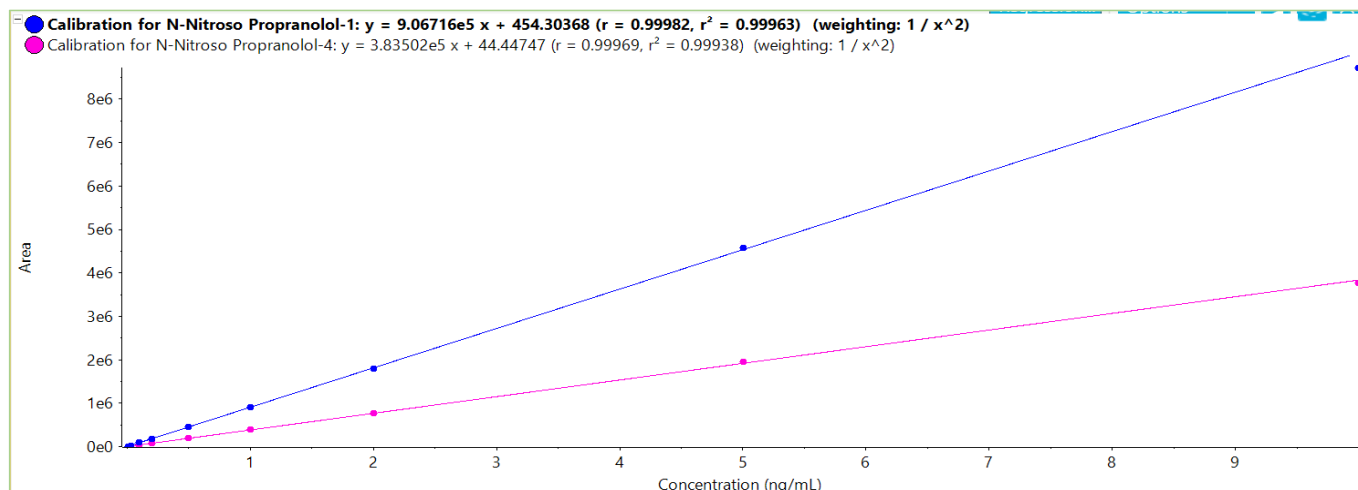


**Figure 3.** Calibration curves were generated using the quantifier (blue) and the qualifier (pink) transition responses. The linear dynamic range (LDR) was 3 orders of magnitude.

Two MRM transitions were used to provide an additional level of specificity for N-nitroso propranolol quantification (Table 3). The ion ratios for the N-nitroso propranolol spike samples are similar to the ratios observed in the standard solution, indicating that under the optimized LC-MS/MS conditions, the interferences at the retention time of N-nitroso propranolol are

negligible. Representative XICs shown in Figure 4 highlight the ion ratio lines with a tolerance of  $\pm 30\%$ .

Recovery was evaluated in 3 matrices, including the placebo, API and drug product. Recovery was assessed at 0.01 ng/mL, 0.03 ng/mL, 1 ng/mL and 5 ng/mL for the placebo and API samples. Since N-nitroso propranolol was detected in the drug product at 0.185 ng/mL, the recovery analysis was performed at 1 ng/mL and 5 ng/mL. In the case of the placebo sample, no area response was observed at the retention time of N-nitroso propranolol. Table 3 summarizes the recoveries observed for all evaluated matrix samples. Recoveries between 85% and 111% were achieved, indicating that the extraction workflow achieved good recovery for N-nitroso propranolol in both the drug substance and product.

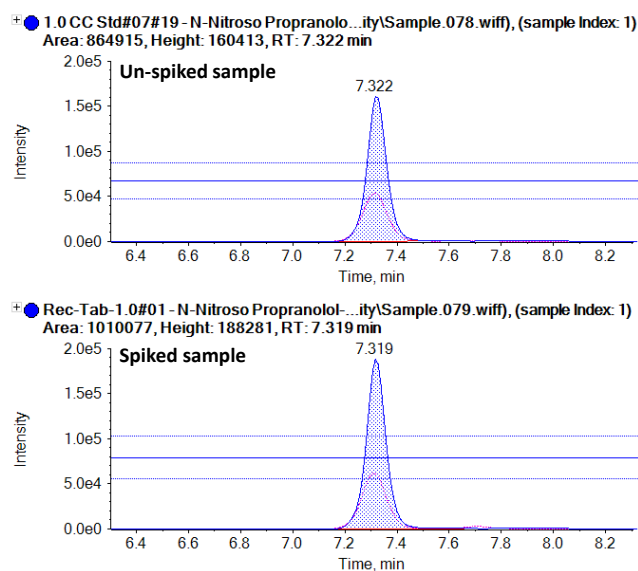

**Figure 4.** XICs of N-nitroso propranolol for standard solutions and spiked samples. Overlaid XICs of N-nitroso propranolol from the standard solution (top) and spiked with N-nitroso propranolol (bottom). The quantifier (blue) and qualifier (pink) ions are overlaid. Ion ratio lines highlight a tolerance of  $\pm 30\%$ , which was set as the optimal ratio between the quantifier and qualifier fragment responses. The standard sample (top) contained 0.185  $\mu\text{g/g}$  of N-nitroso propranolol.

**Table 3. Recovery results from the standard solutions and spiked samples.**

Sample type	Concentration (ng/mL)	Mean area response in control samples (Mean ± SD)	Mean area response in standard solutions (Mean ± SD)	Mean area response in spiked samples (Mean ± SD)	Recovery (%)
Placebo	0.01	Not detected	9979 ± 642	10622 ± 1155	106
	0.03		28446 ± 924	29886 ± 886	105
	1		1005374 ± 36034	934359 ± 6170	92.9
	5		4878501 ± 221584	4858043 ± 37729	99.6
API	0.01	11248 ± 698.7	9979 ± 642	19525 ± 863	91.9
	0.03		28446 ± 924	36689 ± 420	92.4
	1		1005374 ± 36034	945876 ± 5629	93.0
	5		4878501 ± 221584	4651926 ± 28735	95.1
Drug product	1	168938 ± 8961.1	1005374 ± 36034	1008797 ± 8748	85.9
	5		4878501 ± 221584	4669906 ± 45067	92.5


## Compliant-ready SCIEX OS software

SCIEX OS software is a closed system and requires records and signatures to be stored electronically, meeting the regulations outlined in 21 CFR Part 11. SCIEX OS software

can open raw data files from any visible storage location within a closed network by using designated processing workstations. Figure 5 illustrates 3 types of controls required for 21 CFR Part 11 compliance. The workflow presented here is fully compliant with these guidelines, as SCIEX provides 1) technical controls over hardware and software configuration, 2) network security and secure operating systems and policies and 3) procedures and user training (Figure 5).

## Conclusions

- Low-level quantification of the N-nitroso propanol impurity was achieved using the QTRAP 6500+ system
- Accurate and highly reproducible quantification of N-nitroso propanol was demonstrated. Linearity was achieved between 0.01 to 10 ng/mL, generating an LDR of 3 orders of magnitude
- The method demonstrated a simple sample preparation workflow with high recovery. The recovery was evaluated in placebo, API and drug product matrices.
- A single platform for streamlined data processing and management with SCIEX OS software was presented. The SCIEX OS software is also a 21 CFR Part 11-compliant platform that can be easily integrated into a regulated environment.

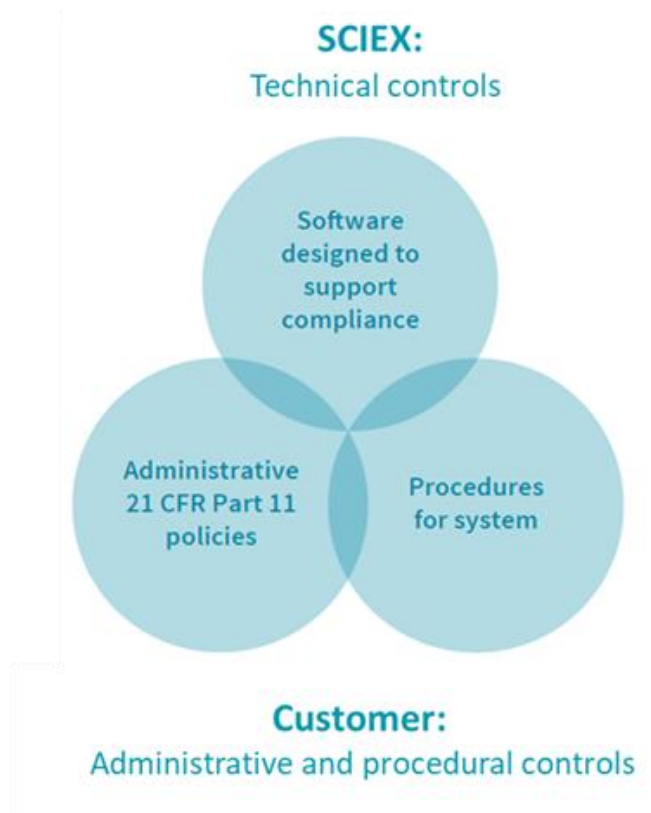**Figure 5. Controls required for 21 CFR Part 11 compliance.**

## References

1. Partani P. et al., Simultaneous determination of propranolol and 4-hydroxy propranolol in human plasma by solid phase extraction and liquid chromatography/electrospray tandem mass spectrometry. [J Pharm Biomed Anal. 2009 Dec 5;50\(5\):966-76.](#)
2. Robbiano L, Martelli A, Allavena A, Mazzei M, Gazzaniga GM, Brambilla G. Formation of the N-nitroso derivatives of six beta-adrenergic-blocking agents and their genotoxic effects in rat and human hepatocytes. [Cancer Res. 1991 May 1;51\(9\):2273-9.](#)
3. Recalls propranolol hydrochloride capsules due to a nitrosamine impurity - Canada.ca  
<https://recalls-rappels.canada.ca/en/alert-recall/pfizer-recalls-inderal-propranolol-hydrochloride-capsules-due-nitrosamine-impurity>
4. Control of nitrosamine impurities in human drugs – guidance for industry, FDA, [February 2021.](#)
5. International Council for Harmonisation of Technical Requirements for Pharmaceuticals for Human Use; “ICH Harmonised Guideline - Assessment And Control Of DNA Reactive (Mutagenic) Impurities in Pharmaceuticals to Limit Potential Carcinogenic Risk, [M7\(R1\)](#)”; [March 31, 2017](#)
6. Brambilla G, Martelli A, Sottofattori E. Nitrosation of propranolol under simulated gastric conditions. Carcinogenesis. [1995 May;16\(5\):1239-42.](#)

The SCIEX clinical diagnostic portfolio is For In Vitro Diagnostic Use. Rx Only. Product(s) not available in all countries. For information on availability, please contact your local sales representative or refer to <https://sciex.com/diagnostics>. All other products are For Research Use Only. Not for use in Diagnostic Procedures.

Trademarks and/or registered trademarks mentioned herein, including associated logos, are the property of AB Sciex Pte. Ltd. or their respective owners in the United States and/or certain other countries (see [www.sciex.com/trademarks](http://www.sciex.com/trademarks)).

© 2023 DH Tech. Dev. Pte. Ltd. MKT-25212-A



### Headquarters

500 Old Connecticut Path | Framingham, MA 01701 USA  
Phone 508-383-7700  
[sciex.com](http://sciex.com)

### International Sales

For our office locations please call the division headquarters or refer to our website at [sciex.com/offices](http://sciex.com/offices)