

Release Notes

LibraryView™ Software

Version 1.0

Features

- The LibraryView software provides workflows for building chemical compound libraries, searching libraries and subsets of libraries, reviewing data to confirm accuracy and quality, and creating reports.
- Enhanced spectral matching (% purity) for MRM (multiple reaction monitoring), EMS (enhanced mass spectrum), and IDA (information dependent acquisition) workflows is also provided.
- Automated searching of multiple .wiff files is supported.

Where to Get Help

Other Documentation

- LibraryView software *Help*
- *LibraryView™ 1.0 Software – Getting Started – Quick Reference Card*
- *LibraryView™ 1.0 Software Installation Guide*

Support

- support@absciex.com
- www.absciex.com
- Customer Documentation: techpubs@absciex.com


Requirements

- Windows XP operating system with SP3 or the Windows 7 (32-bit) operating system, with or without SP1.
- Microsoft Office 2007 or 2010.
- To create .pdf versions of reports when Microsoft Office 2007 is installed, the 2007 Microsoft Add-in: Microsoft Save as PDF must be installed. The free download is available from www.microsoft.com/download.
- To create acquisition and quantitation methods using the LibraryView software, Analyst® software version 1.6 or higher must be installed.

Known Issues and Limitations

- When the LibraryView software is installed on a computer configured with the Windows 7 operating system, the software must be installed by a user with Administrator privileges.

- When the TIC or XIC view is maximized, the zooming feature becomes unresponsive. To unlock the zoom view, click any sample name shown in the Search pane.
- Occasionally, when the user alternates between the TIC and the XIC and then attempts to switch between Experiment 1 and Experiment 2 in the TIC pane, the view becomes locked and does not switch between the chromatograms in Experiment 1 and Experiment 2. To unlock the view, click the Sample Name in the Search pane.
- Occasionally, the TIC does not display correctly (displays spectral-like data rather than chromatographic peaks). To display the TIC correctly, click anywhere inside the Search pane to refresh the chromatogram.
- When the Acquired Spectrum pane is maximized, the zooming feature becomes unresponsive. To unlock the zoom view, click anywhere inside the Searched Spectra pane.
- Users may find that searching spectra is slower when the Grid and Stack views are activated.
- When performing a library search on .wiff files generated using the AB SCIEX TripleTOF[®] 5600+ system, the software defaults to a Collision Energy of zero regardless of the value set in the Standard Parameters section of the New Search dialog. To make sure that search results are obtained, clear the Collision Energy Tolerance check box before beginning the search.

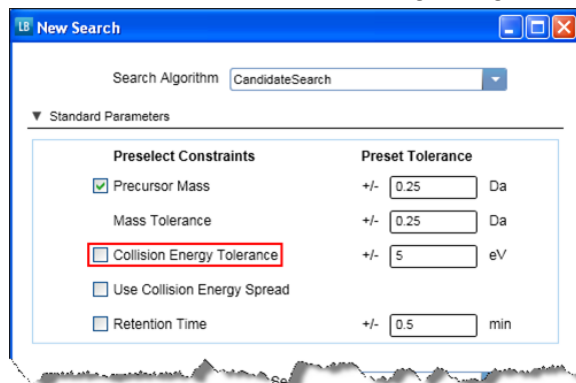

- Refer to the *LibraryView™ 1.0 Software Installation Guide* for a list of potential software installation issues and the recommended solutions.

This document is provided to customers who have purchased AB SCIEX equipment to use in the operation of such AB SCIEX equipment. This document is copyright protected and any reproduction of this document or any part of this document is strictly prohibited, except as AB SCIEX may authorize in writing.

Software that may be described in this document is furnished under a license agreement. It is against the law to copy, modify, or distribute the software on any medium, except as specifically allowed in the license agreement. Furthermore, the license agreement may prohibit the software from being disassembled, reverse engineered, or decompiled for any purpose.

Portions of this document may make reference to other manufacturers and/or their products, which may contain parts whose names are registered as trademarks and/or function as trademarks of their respective owners. Any such usage is intended only to designate those manufacturers' products as supplied by AB SCIEX for incorporation into its equipment and does not imply any right and/or license to use or permit others to use such manufacturers' and/or their product names as trademarks.

AB SCIEX makes no warranties or representations as to the fitness of this equipment for any particular purpose and assumes no responsibility or contingent liability, including indirect or consequential damages, for any use to which the purchaser may put the equipment described herein, or for any adverse circumstances arising therefrom.

For Research Use Only. Not for use in diagnostic procedures.

The trademarks mentioned herein are the property of AB Sciex Pte. Ltd. or their respective owners.

AB SCIEX™ is being used under license.

AB SCIEX

71 Four Valley Dr., Concord, Ontario, Canada. L4K 4V8.

AB SCIEX is ISO 9001 registered.

© 2012 AB SCIEX.

Printed in Canada.

