




# High-resolution, specific analysis of arginine and citrulline on the Echo<sup>®</sup> MS+ system with ZenoTOF 7600 system

Aaron Stella, Jacob W. McCabe and Anuja Bhalkikar  
SCIEX, USA

This technical note demonstrates a high-throughput, high-resolution method for differentiating amino acids that differ by only 1 Dalton. Distinguishing 2 compounds that differ by 1 Dalton can be challenging on a low-resolution instrument, as mass interference is likely to occur and the quantitative values obtained from these instruments will not be trustworthy. Arginine and citrulline are amino acids that differ in molecular weight by 1 Dalton. Differentiating between these amino acids is critical for identifying them as biomarkers or enzyme activators. Chromatographic methods used to separate these compounds can be complex and time consuming. In this work, the Echo<sup>®</sup> MS+ system with ZenoTOF 7600 system was used to resolve arginine from citrulline in a high-throughput manner at a rate of 5 seconds per sample. A Zeno MRM<sup>HR</sup> scan was employed to generate high-resolution product ions for arginine and citrulline simultaneously (Figure 1). Data review was conducted using the SCIEX OS software with library searching for added confidence in identifying arginine and citrulline mixed in solution.

## Key features of the rapid, high-resolution analysis of arginine and citrulline

- **Accurate and specific quantitation of arginine (175.1 Da/60.0562 Da) and citrulline (176.1 Da/159.0960 Da):** Quantify compounds that are 1 Dalton apart with unique product ions and eliminate issues from mass interference
- **Rapid analysis:** Process Zeno MRM<sup>HR</sup> data performed at 5 seconds per sample to achieve ample scans for quantitation
- **Linear and accurate calibration curves:** Quantify values with a wide dynamic range, spanning 3 orders of magnitude from 1.95 μM to 1000 μM with r values >0.99, and accuracies within ±20% of their stated concentrations
- **Full scan and targeted MS/MS data:** TOF and Zeno MRM<sup>HR</sup> data can be collected within the same full scan MS experiment with targeted MS/MS of precursor masses

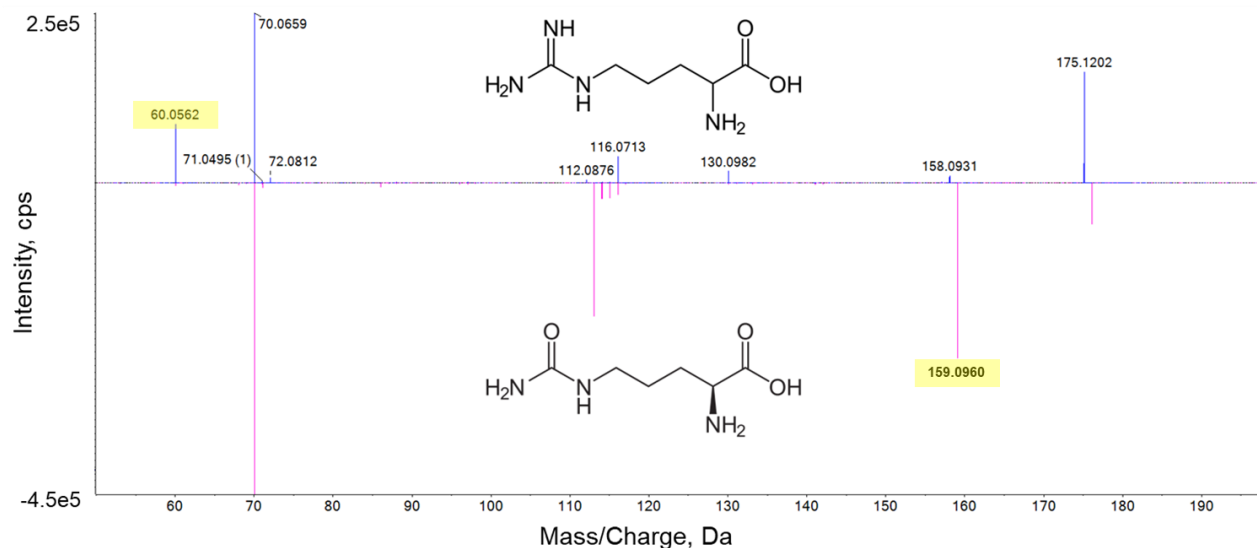

**Figure 1. Arginine and citrulline product ions.** The product ions for arginine (blue trace) and citrulline (pink trace, inverted) were simultaneously obtained from a 5-second Zeno MRM<sup>HR</sup> scan of a sample containing arginine and citrulline. The yellow highlighting indicates the unique product ions that occurred at 60.0562 Da for arginine and 159.0960 Da for citrulline.

## Introduction

Arginine and citrulline are amino acids which only differ in molecular weight by 1 Dalton. Many analytical methods require these compounds to be derivatized for their analysis due to the lack of a natural chromophore in their structure.<sup>1</sup> Derivatization adds time to sample preparation and subsequent liquid chromatography methods can take minutes to complete.<sup>2</sup> These amino acids can be used as biomarkers and enzyme activators but their proximity in molecular weight makes them difficult to be accurately quantified even if derivatization is performed. Here, we present a rapid, high-resolution solution to identify and quantify arginine and citrulline.

## Methods

**Sample preparation:** Arginine, citrulline and a combination of equimolar arginine and citrulline were serially diluted in water to a concentration range of 1.95 $\mu$ M to 1000 $\mu$ M. All samples contained 250 $\mu$ M of phenylalanine as an internal standard (IS).

**Acoustic ejection:** A total of 70 nL of sample was ejected in 5 second intervals at 10 Hz in wide peak mode. Methanol with 0.1% formic acid was the carrier solvent and a flow rate of 400  $\mu$ L/min was used.

**Mass spectrometry:** A Zeno MRM<sup>HR</sup> method was used to quantify arginine and citrulline using both a common product ion and a unique product ion. The data from the common and unique product ions were then compared and the amount of mass interference was calculated. A single unique product ion for the IS was analyzed (Tables 1-4).

**Data processing:** The Analytics and Explorer modules of SCIEX OS software were used for data processing. The NIST library was used to ensure that the analyte identification was correct for arginine and citrulline. The amount of mass interference was

expressed in percent and was calculated using the following equation:

$$\text{Mass interference (\%)} = \left(\frac{A}{B}\right) \times 100$$

where A = the calculated concentration of the interfering analyte and B = the actual concentration of the targeted analyte. Both concentrations were expressed in  $\mu$ M.

**Table 1. Source parameters and values.**

| Parameter                   | Value    |
|-----------------------------|----------|
| Polarity                    | Positive |
| Spray voltage (V)           | 5500     |
| Curtain gas (psi)           | 35       |
| CAD gas (psi)               | 11       |
| Ion source gas 1 (psi)      | 90       |
| Ion source gas 2 (psi)      | 75       |
| Temperature ( $^{\circ}$ C) | 400      |

**Table 2. TOF MS parameters and values.**

| Parameter                  | Value      |
|----------------------------|------------|
| Scan type                  | Zeno MRMHR |
| TOFMS start mass (m/z)     | 50         |
| TOFMS stop mass (m/z)      | 200        |
| Accumulation time (s)      | 0.1        |
| Declustering potential (V) | 60         |
| Time bins to sum           | 4          |

**Table 3. TOF MS/MS parameters and values.**

| Parameter            | Value |
|----------------------|-------|
| Q1 resolution        | Unit  |
| Zeno pulsing         | On    |
| Zeno threshold (cps) | 2000  |

**Table 4. Zeno MRM<sup>HR</sup> parameters and values.**

| Compound ID        | Precursor ion | TOF start (m/z) | TOF stop (m/z) | Accumulation time (s) | Declustering potential (V) | Collision energy (V) | CE spread (V) | Time bins to sum |
|--------------------|---------------|-----------------|----------------|-----------------------|----------------------------|----------------------|---------------|------------------|
| Arginine           | 175.1         | 50              | 200            | 0.07                  | 50                         | 23                   | 10            | 4                |
| Citrulline         | 176.1         | 50              | 200            | 0.07                  | 50                         | 23                   | 10            | 4                |
| Phenylalanine (IS) | 166.2         | 50              | 200            | 0.07                  | 50                         | 21                   | 10            | 4                |

## High-resolution product ions

The Zeno MRM<sup>HR</sup> scan provides a TOF MS scan in addition to performing an MRM scan for precursor masses. The user enters a precursor ion for each analyte or IS to be measured and then the user can choose the TOF MS/MS mass range to be scanned. The collision energy can be set to a fixed number or spread around a central collision energy value.

When combined with wide peak mode ejections, the resulting scans include data-rich product ions that were individually selected for analyte and IS processing with a library search using the Analytics module in SCIEX OS software (Figure 1). Wide peak mode with a 5 second sample ejection interval was chosen to ensure that 12 scans were obtained for each sample (Figure 2).

Table 5. Zeno MRM<sup>HR</sup> transitions selected for quantitation.

| Analyte    | Common product ion (m/z) | Unique product ion (m/z) |
|------------|--------------------------|--------------------------|
| Arginine   | 175.1256/70.0659         | 175.1256/60.0562         |
| Citrulline | 176.1094/70.0659         | 176.1094/159.0960        |

## Calibration curves

Calibration curves with common and unique product ions were constructed from the samples containing a combination of arginine and citrulline at concentrations ranging from 1.95 $\mu$ M to 1000 $\mu$ M (Figure 3).


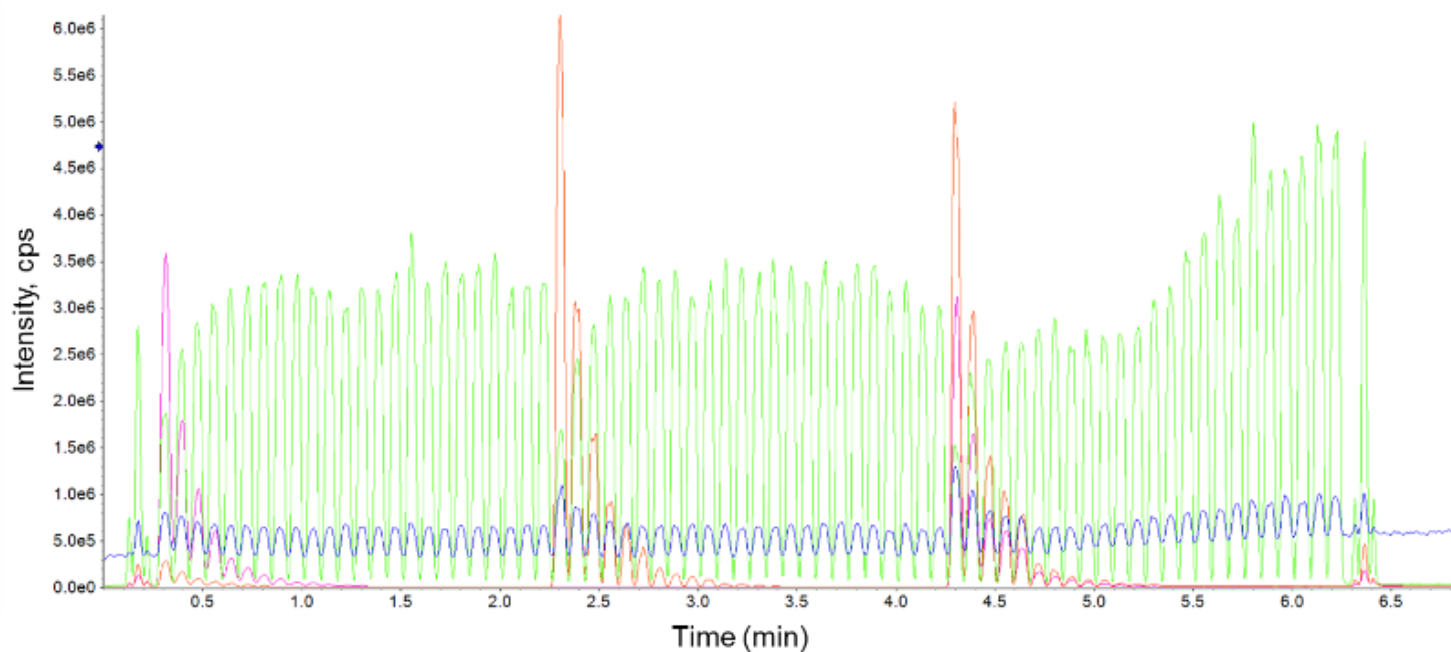


Figure 2. Rapid, quantitative analysis of arginine, citrulline and phenylalanine using the Echo<sup>®</sup> MS+ system with ZenoTOF 7600 system. The sample was analyzed in triplicate. Arginine is shown in pink and citrulline is shown in orange. Phenylalanine was used as an IS and is shown in green. The TOF MS scan is shown in blue.

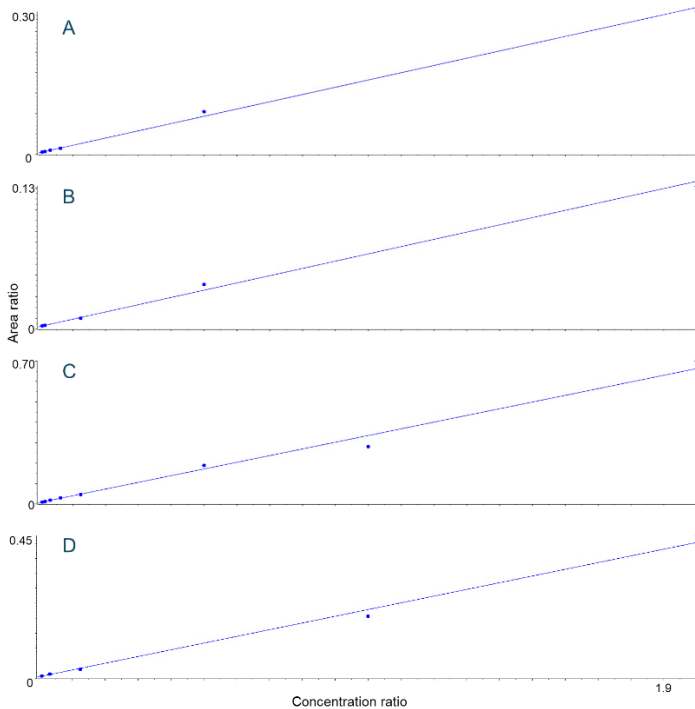

**Figure 3. Calibration curves for the analysis of arginine and citrulline.** Calibration curves for arginine based on a common (A) or unique product ion (B) are shown, in addition to calibration curves for citrulline based on a common (C) or unique product ion (D). The area ratio was calculated by dividing the analyte peak area by the IS peak area. The concentration ratio was calculated by dividing the analyte concentration by the IS concentration.

The calibration curves were then used to quantify arginine and citrulline and were applied to calculate the amount of mass interference. All calibration point accuracies were within  $\pm 20\%$  of their assigned values and all calibration curve  $r$  values were  $>0.994$ .


### Assessment of mass error

When quantifying with unique product ions for arginine and citrulline, no mass interference was detected for either analyte. A small percentage of mass interference was detected when a common product ion was used for quantitation. Only in high concentration arginine samples was arginine mis-identified as citrulline. Less than 2% mass interference was detected in the citrulline common product ion channel when testing a 1000 $\mu\text{M}$  and a 500 $\mu\text{M}$  arginine sample (Table 6).

**Table 6. Percent (%) mass interference detected in high concentration samples of arginine and citrulline.**


| Analyte, Concentration         | Arginine, Common Product Ion | Arginine, Unique Product Ion | Citrulline, Common Product Ion | Citrulline, Unique Product Ion |
|--------------------------------|------------------------------|------------------------------|--------------------------------|--------------------------------|
| Arginine, 1000 $\mu\text{M}$   | <i>n/a</i>                   | <i>n/a</i>                   | 1.72                           | 0                              |
| Citrulline, 1000 $\mu\text{M}$ | 0                            | 0                            | <i>n/a</i>                     | <i>n/a</i>                     |
| Arginine, 500 $\mu\text{M}$    | <i>n/a</i>                   | <i>n/a</i>                   | 0.407                          | 0                              |
| Citrulline, 500 $\mu\text{M}$  | 0                            | 0                            | <i>n/a</i>                     | <i>n/a</i>                     |

No mass interference was detected in samples with concentrations  $\leq 250\mu\text{M}$  of either analyte, independent of the type of product ion used for quantitation. With regards to the small amount of mass interference observed in the arginine sample when quantifying with the common product ion for citrulline, the library hit yielded “no match,” resulting in a library score of 0 (Figure 4).

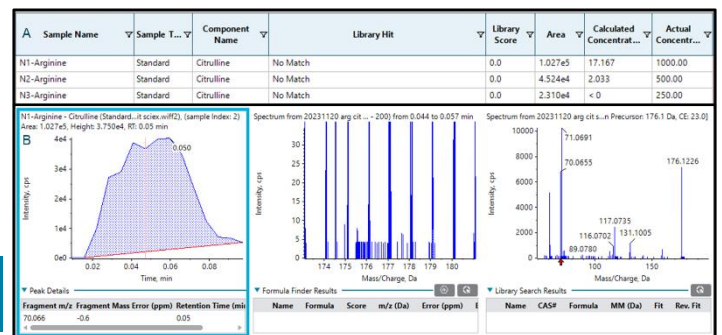

**Figure 4. Samples containing arginine at 1000 $\mu\text{M}$ , 500 $\mu\text{M}$  and 250 $\mu\text{M}$  analyzed using the citrulline common product ion MRM channel.** This analysis registered 17.17 $\mu\text{M}$ , 2.033 $\mu\text{M}$  and  $< 0\mu\text{M}$  of mass interference, respectively, indicated in the “calculated concentration” column (A). The extracted ion chromatogram peak, TOF MS and TOF MS/MS scan with library search results are shown for the 1000 $\mu\text{M}$  arginine sample (B). When quantified with the citrulline common product ion, the results showed 17.17 $\mu\text{M}$  mass interference.

The peak area ratios obtained from the citrulline common product ion MRM channel were plotted against the peak area ratios obtained from the citrulline unique product ion MRM channel in samples with arginine concentrations of 1000 $\mu$ M, 500 $\mu$ M and 250 $\mu$ M (Figure 5). The mass interference observed with the common product ion was neutralized when the same 3 samples were analyzed using the citrulline unique product ion.

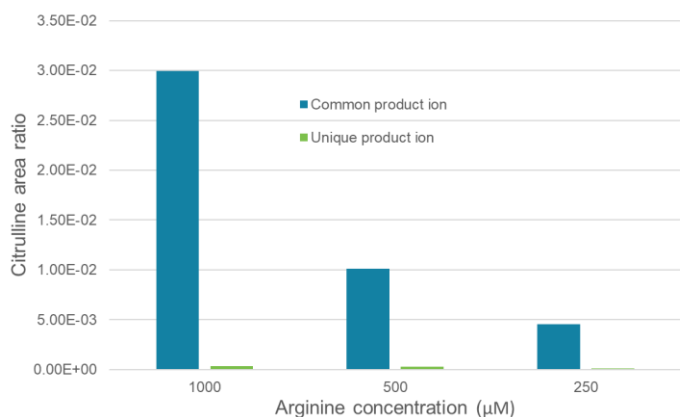

**Figure 5. Observed arginine peak area ratios determined by the citrulline common and unique product ions.** The peak area ratios were determined at concentrations of 1000 $\mu$ M, 500 $\mu$ M and 250 $\mu$ M.

Quantitation of arginine and citrulline was accurate and specific based on the calibration curve data. To further ensure accurate identification of arginine and citrulline, we utilized the NIST library. At concentrations as low as 1.95 $\mu$ M, the NIST library accurately identified the correct compounds in the samples (Figure 6). For the sample containing arginine, only arginine and the phenylalanine IS were detected and for the sample containing citrulline, only citrulline and the phenylalanine IS were detected. For the combined arginine and citrulline sample, arginine, citrulline and the phenylalanine IS were all detected. All library scores suggested a good level of confidence, as the lowest library score was 79.9.

| Sample Name                 | Component Name   | Fragment Mass | Library Hit                                  | Library Score |
|-----------------------------|------------------|---------------|--|---------------|
| N10-Arginine                | Arginine         | 70.066        | DL-Arginine (NIST) [Smart Confirmation]      | 84.9          |
| N10-Arginine                | Citrulline       | 70.066        | No Match                                     | 0.0           |
| N10-Arginine                | Phenylalanine IS | 120.082       | DL-Phenylalanine (NIST) [Smart Confirmation] | 100.0         |
| N10-Arginine                | Arginine 2       | 60.056        | DL-Arginine (NIST) [Smart Confirmation]      | 79.9          |
| N10-Arginine                | Citrulline 2     | 159.096       |  | N/A           |
| O10-Citrulline              | Arginine         | 70.066        |  | N/A           |
| O10-Citrulline              | Citrulline       | 70.066        | L-Citrulline (NIST) [Smart Confirmation]     | 98.4          |
| O10-Citrulline              | Phenylalanine IS | 120.082       | DL-Phenylalanine (NIST) [Smart Confirmation] | 100.0         |
| O10-Citrulline              | Arginine 2       | 60.056        |  | N/A           |
| O10-Citrulline              | Citrulline 2     | 159.096       | L-Citrulline (NIST) [Smart Confirmation]     | 98.9          |
| P10-Arginine and Citrulline | Arginine         | 70.066        | L-Arginine (NIST) [Smart Confirmation]       | 88.5          |
| P10-Arginine and Citrulline | Citrulline       | 70.066        | L-Citrulline (NIST) [Smart Confirmation]     | 96.5          |
| P10-Arginine and Citrulline | Phenylalanine IS | 120.082       | DL-Phenylalanine (NIST) [Smart Confirmation] | 100.0         |
| P10-Arginine and Citrulline | Arginine 2       | 60.056        | L-Arginine (NIST) [Smart Confirmation]       | 84.3          |
| P10-Arginine and Citrulline | Citrulline 2     | 159.096       | L-Citrulline (NIST) [Smart Confirmation]     | 97.6          |

**Figure 6. Library search confirmation performed using the Analytics module in SCIEX OS software.** The sample name, MRM component, fragment mass, library hit and library score for samples containing arginine, citrulline and a 1.95 $\mu$ M solution containing arginine and citrulline are shown. Each sample included 250 $\mu$ M phenylalanine as an IS.

Mass interference when quantifying with the common product ion for citrulline was only observed at concentrations >250 $\mu$ M. For quantitation at concentrations <250 $\mu$ M, either the unique or common MRM channel could produce specific results. Ultimately, choosing a unique product ion is preferred, as the likelihood of registering mass interference is much lower, given the product ion is exclusive to the analyte of interest. However, these results show that when unique product ions are not available, the common product ions from the high-resolution, Echo<sup>®</sup> MS+ system with ZenoTOF 7600 system could be used for sensitive and specific quantitation if excessively high concentrations of potential mass interfering substances are avoided.

## Conclusions

- Achieved specific quantitation of arginine and citrulline, which are 1 Dalton apart in parent mass
- Wide peak mode provided 10 data points per peak with a 3.7 ms ejection of 70 nL of sample
- Rapid, linear quantitation with an IS across 3 orders of magnitude
- Mass interference was significantly reduced by quantifying with a unique product ion in Zenon MRM<sup>HR</sup> mode
- Mass interference was absent when quantifying with unique product ions
- NIST library search was used to increase our identification confidence

## References

1. H. Mao, W. Wei, W. Xiong, Y. Lu, B. Chen, Z. Liu. Simultaneous Determination of L-Citrulline and L-Arginine in Plasma by High Performance Liquid Chromatography. *Clinical Biochemistry*. September 2010, 43(13-14);1141-7.  
<https://www.sciencedirect.com/science/article/abs/pii/S0009912010002687?via%3Dihub>
2. Perez-Neri, E. Castro, S. Montes, M-C. Boll, J. Barges-Coll, JL, Soto-Hernandez, C. Rios. Arginine, Citrulline and Nitrate Concentrations in the Cerebrospinal Fluid from Patients with Acute Hydrocephalus. *Journal of Chromatography B*. 2007, 851;250-6.  
<https://www.sciencedirect.com/science/article/abs/pii/S1570023206008658?via%3Dihub>

The SCIEX clinical diagnostic portfolio is For In Vitro Diagnostic Use. Rx Only. Product(s) not available in all countries. For information on availability, please contact your local sales representative or refer to <https://sciex.com/diagnostics>. All other products are For Research Use Only. Not for use in Diagnostic Procedures.

Trademarks and/or registered trademarks mentioned herein, including associated logos, are the property of AB Sciex Pte. Ltd. or their respective owners in the United States and/or certain other countries (see [www.sciex.com/trademarks](http://www.sciex.com/trademarks)).

Echo<sup>®</sup> and Echo<sup>®</sup> MS are trademarks or registered trademarks of Labcyte, Inc. in the United States and other countries, and are being used under license.

© 2024 DH Tech. Dev. Pte. Ltd. MKT-27824-A



### Headquarters

500 Old Connecticut Path | Framingham, MA 01701  
USA  
Phone 508-383-7700  
[sciex.com](http://sciex.com)

### International Sales

For our office locations please call the division headquarters or refer to our website at [sciex.com/offices](http://sciex.com/offices)