

Mass Spec Made Simple

The SCIEX Topaz™ System for Clinical Diagnostics

Harness the power of mass spectrometry...while leaving behind the complexity


The intrinsic complexity of traditional LC-MS/MS systems has historically made this powerful technology inaccessible to many busy clinical labs with constrained resources. The result has been an undesirable compromise of costly send-out testing, extended turn-around times, and limiting assay performance. The new SCIEX Topaz System now bridges this gap, with novel ClearCore™ MD Software. The intuitive design of ClearCore MD Software simplifies and streamlines workflows and method development, bringing the power of SCIEX LC-MS/MS technology to the hands of the routine user.

Everything at your fingertips

- Intuitive, tile-based user interface
- Access to all system modules with a single click

Each module offers

- Simplified menus
- Select number of options
- Streamlined navigation for common tasks

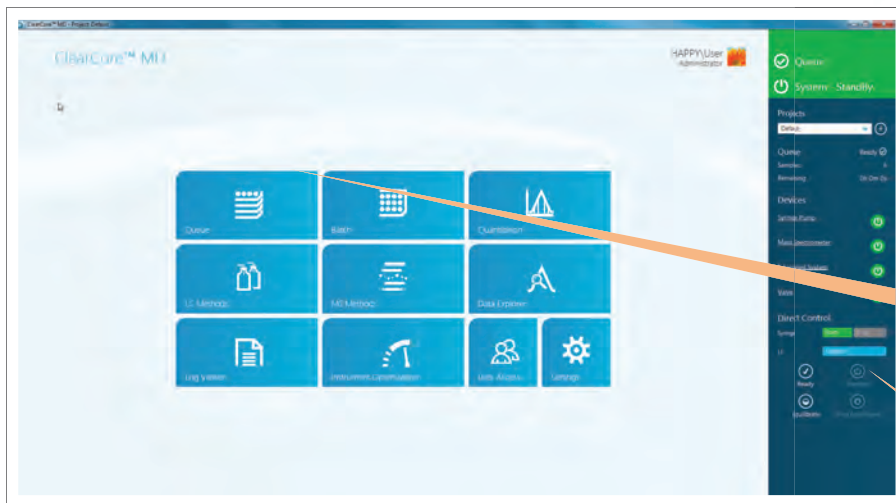

Access to complete functionality via a simple tile-based grid

Instant view of sample and instrument status

All systems go

An Instrument/Sample Status Panel can be viewed from any ClearCore™ MD module the user is working. Key features of this status panel include:

- **Visual indicator lights to instantly inform the user of current system state**
- **Instant status updates from anywhere in the software, so no time is wasted**

Take control

ClearCore MD Software offers direct device control via the status panel or individual modules. Features include:

- **Component specific tabs to simplify control**
- **Visual tools and real time feedback**
- **Rapid method development and deployment by user**
- **An end to complex control panel navigation**

Method development simplified


ClearCore MD Software provides enhanced ease of use for rapid method development, including

- **Flexibility of discrete LC and MS method modules**
- **Automatic Optimization with minimal user input**
- **Automatic optimization report generation.**

Stay connected and secure


ClearCore MD Software provides many features to maintain integrity of sample analysis

- **Multi-level user access**
 - **Functionality determined by user level**
- **Bidirectional LIS communication**
 - **Transcription-free import and export of sample information and results**

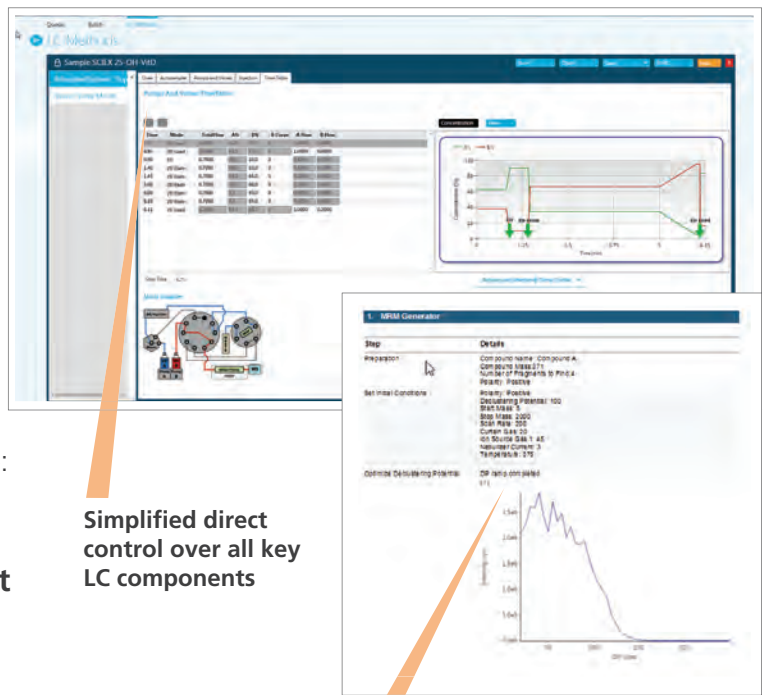

Simplified direct control over all key LC components

Automated MRM optimization for rapid method development by user

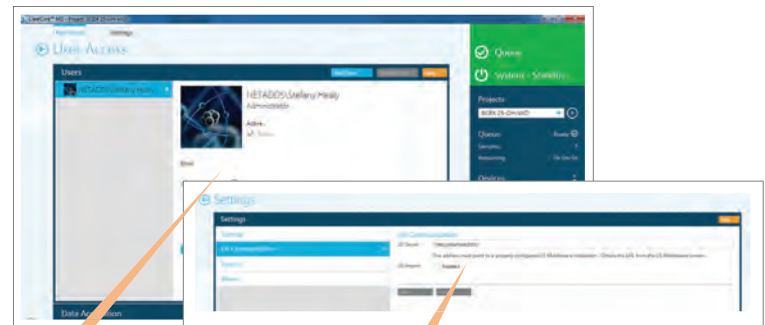

Simple restriction of user access via pre-defined user privilege levels

Easy setup of LIS communication for sample batch setup and results export



Learn more about the Topaz system at sciex.com/MassSpecMadeSimple